



ANDERSEN
STERILIZERS

AN74 Series Features & Specifications



**Protect your patients,
Preserve your instruments**

The most cost-effective, gentle and
compatible room-temperature
tabletop sterilizer.

Uses **90% less EO**
than any other system.

Anprolene AN74i

Suggested use:

Space-saving design with
small practices in mind.
Knockout your sterilization
needs, not your budget.

Chamber dimensions:

19.25"Lx13.75"Wx12.75"H

Anprolene AN74j

Suggested use:

Often found in zoos and
equine practices due to
its ability to sterilize
oversized instruments
and large loads.

Chamber dimensions:

29.25"Lx17.75"Wx12.5"H

Anprolene AN74ix

Suggested use:

Suitable for general use
and medium-volume
veterinary clinics. Can
accommodate endoscopes
and long instruments.

Chamber dimensions:

29.25"Lx13.75"Wx12.75"H



17.6 GRAM AMPOULE



35 LITER LINER BAG



FREE OPERATOR TRAINING

Bring the proven reliability of EO sterilization into
your facility with the Anprolene System.

The Anprolene system has always offered the **most
gentle 100% EO sterilization cycle** on the market:
Room temperature sterilization (minimum 68° F),
no steam injection, no deep vacuum, and **only 17.6
grams of gas**.

Unmatched sterilization compatibility. Protects
your most **delicate instruments, electronics, drills,
fiber optics, hand pieces and cameras**. As well as
those **plastic, fabric, cellulose or rubber items**
that you can't sterilize with other methods.

FREE operator training for the life of your sterilizer.



MOST COMPATIBLE

STERILIZER IS INCLUDED WITH A STARTER KIT :

Starter Kit - All you need for 28 cycles!

- 2 cases 17.6g Ampoules and Sterilization Liner Bags / **AN7916**
- 2 pouches Humidichip RH Stabilizer Chips / **AN1071**
- 1 unit Humiditube Holder / **AN1072**
- 2 boxes Dosimeter Chemical Indicators / **AN87**
- 1 box EO Indicator Tape / **AN2250**
- 1 box EO Exposure Indicator Strips / **AN85**
- 1 box EZTest Biological Indicators / **AN2203**
- 1 kit Kem Medical Vapor-Trak Monitors / **AN94**
- 1 unit Installation Kit / **AN7401**

Affordable: Our Anprolene line is quite simply the most affordable on the market – both in terms of cost to purchase and per cycle.

Easy Installation: Requires only a 120v or 240v outlet and a 1" exhaust line.



Sterilization and Aeration in the same cabinet: Andersen's purge probe flushes the liner bag with a constant flow of fresh air at the end of the cycle – no need to transfer items to a separate area. Two hours of aeration standard. EO-absorbent items may require additional aeration. Active aeration continues indefinitely until products are removed.

Proven Reliability: Anprolene is used in over 5,000 facilities worldwide. Trusted for more than 60 years. Find out how an Anprolene sterilizer can improve your practice today.

Make Your System Zero Emissions: Andersen abators are a simple, cartridge-based system that employ a dry catalyst resin. The resin converts ethylene oxide to biodegradable organic compounds. Replaceable abator cartridge has a 200 cycle capacity. Spent cartridges are non-hazardous.



ABATOR WITH CARTRIDGE

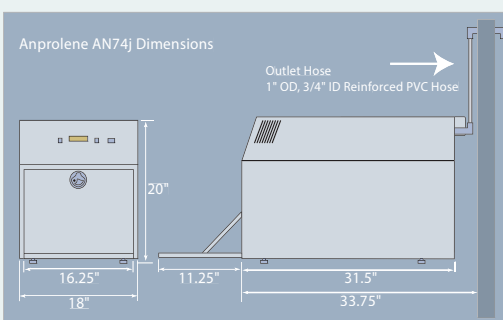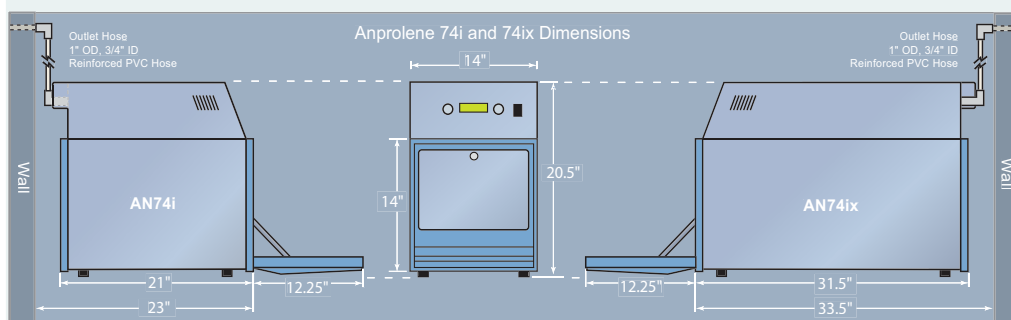

Specifications: AN74i

Model: AN74.64
Interior dimensions: 19.25"L x 13.75"W x 12.75"H
Exterior dimensions: 21"L x 14"W x 20.5"H
Cabinet capacity: 55 l
Total weight: 50.5 lbs.
Power Supply: 115v +/- 10%, 60 Hz/150 Watts
230v available (AN74.65)
Installation Category II

Specifications: AN74ix

Model: AN74.84
Interior dimensions: 29.25"L x 13.75"W x 12.75"H
Exterior dimensions: 31.5"L x 14"W x 20.5"H
Cabinet capacity: 84 l
Total weight: 60.5 lbs.
Power Supply: 115v +/- 10%, 60 Hz/150 Watts
230v available (AN74.85)
Installation Category II

Specifications: AN74j

Model: AN74.74
Interior dimensions: 29.25"L x 17.75"W x 12.5"H
Exterior dimensions: 31.5"L x 18"W x 20"H
Cabinet capacity: 110 l
Total weight: 85 lbs.
Power supply: 115 +/- 10%, 50/60 Hz/150 Watts
230v available (AN74.75)
Installation Category II

PN4428_033121